

SYSTEM 3

Stretch-Glide Mattress

Exercise excellent pressure reduction by additional
Two types property urethane foam

Guidelines on Probable Users (Risk of bedsores)			
None	Low	Medium	High
Stretch-Slim			
Stretch-Fit			
Stretch-Glide			

Upper layer



Top-level pressure reduction for stationary mattresses*.

[Two layers urethane foam with unique cutting treatment]

*Compare with Paramount Bed's mattress lineup

New mechanism of "Glide sheet processing"


Employs "Glide sheet processing" to treat urethane foam and apply a smooth. It reduces a risk of pressure ulcer because of keeping the person from sliding when back raising.

Side edges work to easily stabilize sitting on side.

Use firmer urethane foam on both side of the mattress makes stable sitting posture on side and getting into or out of the bed.

Exercise excellent pressure reduction by stretch system and Urethane foam


Exercise top-level pressure reduction for stationary mattresses by urethane foam and cuts introduced to support certain parts of the body.

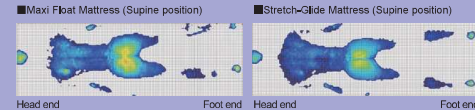


Low Pressure High
(H)175cm, (Weight)65kg, Male
*according to Paramount Bed standards and procedures.




A Upper body part cutting
Correspond spinal column and shoulders when lateral recumbent position.

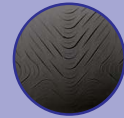


B Buttocks part cutting
Correspond sacrum and ilium, greater trochanter when lateral recumbent position.

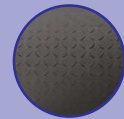


C Legs part cutting
Correspond shape of heels

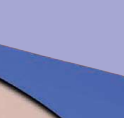

D Heel part
Disperse pressure for heel to calf by smooth slope on lower leg part.

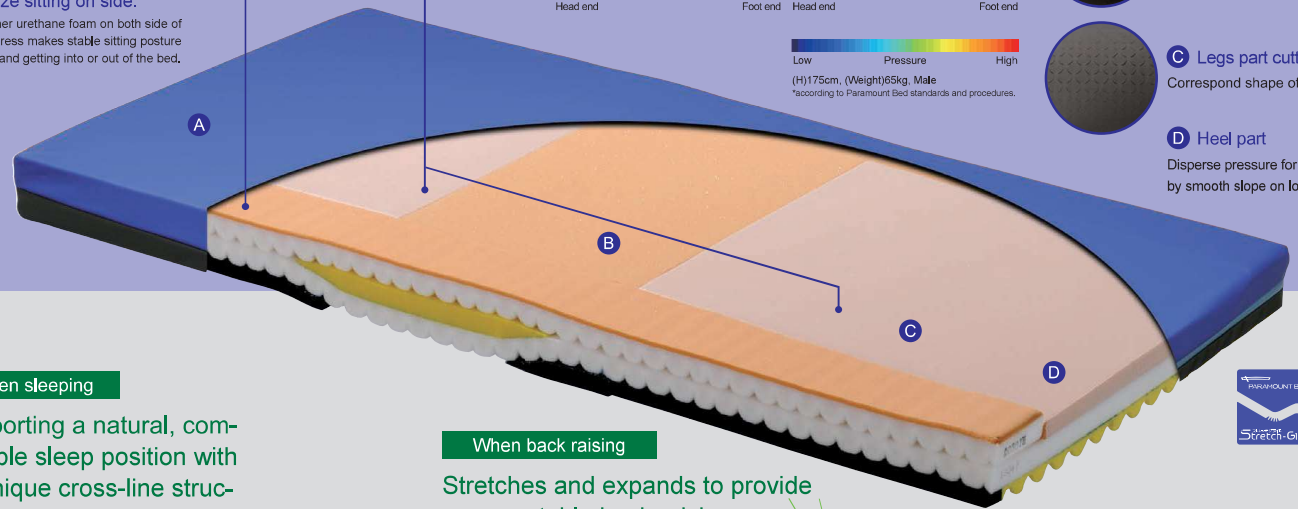

Lower layer

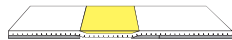


Mattress stretches to reduce patient migration.

[New mechanism: Stretch mattress system]

When sleeping


Supporting a natural, comfortable sleep position with its unique cross-line structure.

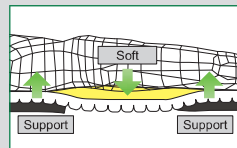

● Providing support for the lumbar, posterior, and femoral regions by area.

When back raising

Stretches and expands to provide a more stable back raising.

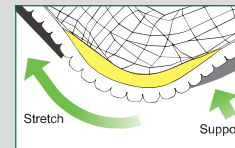

● Stretch with Wave-patterned structure.
● Anchor support keeps posture.



Designed to reduce patient migration

Together with Smart Raise function of bed, it reduces shear force and feeling of oppression during back raising by stretch structure.

Water proof type

Correspond to "Sodium hypochlorite" which disinfects against a wide range of types of bacteria and viruses.

Size	Model No.	Size/Product weight
83cm width	KE-793S	(W)83 x (L)191 x (TH)12.5cm, (Weight)8.5kg
91cm width	KE-791SWA	(W)91 x (L)191 x (TH)12.5cm, (Weight)9kg
86cm width	KE-798SYA	(W)86 x (L)200 x (TH)12.5cm, (Weight)9kg



Washable type

Can be washed with cover

Size	Model No.	Size/Product weight
83cm width	KE-793T	(W)83 x (L)191 x (TH)12.5cm, (Weight)8kg
91cm width	KE-791T	(W)91 x (L)191 x (TH)12.5cm, (Weight)9kg

